




**CAUTION:** This temperature sensor is sensitive to electrostatic discharge (ESD). Use ESD precautionary procedures when handling, or making mechanical or electrical connections to this device in order to avoid performance degradation or loss of functionality.

Three aspects of using a cryogenic temperature sensor are critical to its optimum performance. The first involves the proper mounting of the sensor package, the second relates the proper joining of sensor lead wires and connecting wires. The final concern is the thermal anchoring of the lead wires. Although the sequence in which these areas should be addressed is not fixed, all elements covered under each aspect should be adhered to for maximum operating capabilities of the sensor.

## SENSOR MOUNTING

1. Mounting area should be prepared and cleaned with a solvent such as acetone followed by an isopropyl alcohol rinse. Allow time for the solvents to evaporate before sensor mounting.
2. The list below provides brief instructions on mounting a sensor using a number of different methods. The constraints of your application should dictate the most appropriate mounting method to follow.

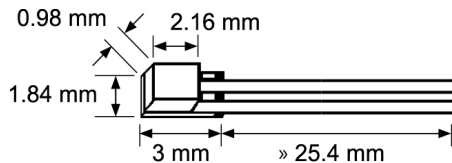

Leads: 32 AWG (0.2 mm diameter)  
Gold-plated Copper (newer models)  
Phosphor-Bronze (older models)  
Sensor Mass: 0.03 gram

**Mechanical** – The preferred method for mechanically mounting a SD package is using the Lake Shore spring loaded clamp. This clamp should be ordered at the time the sensor is ordered (CX-10XX-CO). The clamp holds the SD package in contact with the surface and also allows the sensor to be easily changed or replaced. A thin layer of Apiezon® N grease ( $\leq 0.055$  mm) or a flat 100% indium preform should be used between the sensor and mounting surface to enhance thermal contact. The spring prevents crushing the sensor.

**Indium Solder (100% In)** – A low wattage heat source should be used and the sensor must never exceed 200 °C. The mounting surface and sensor should be tinned with a rosin flux (type RMA is recommended) prior to mounting the sensor. A thin, uniform layer of indium solder should be the goal. Clean both the sensor and mounting surface of residual flux using rosin residue remover. Once the surface area is dry, reheat the mounting surface to the melting point of the solder (156 °C). Press the sensor into position and allow it to warm to the melting point of the solder. Remove heat source and allow sufficient time for the solder to solidify (typically 2 to 3 s) before removing it.

**Apiezon® N Grease** – Is best used as a thermal conductor when the sensor is mounted in a hole or recess and when the sensor is intended to be removed. The sensor should be surrounded with thermal grease and placed into the mounting position. When the temperature is lowered the thermal grease will harden, giving good support and thermal contact.

**IMI 7031 Varnish** – Prepare varnish and apply a thin layer on the mounting surface. Press the sensor firmly against the varnish during curing to ensure a thin bond layer and good thermal contact. Varnish will air dry in 5 to 10 min. Sufficient time must be allowed for the solvents in the varnish to evaporate. There is a small probability of ionic shunting across the sensor during the full cure period of the varnish (typically 12 to 24 h).

**Stycast® 2850FT Epoxy** – Prepare epoxy and apply a thin layer on the mounting surface. Press the sensor firmly into the epoxy during curing to assure a thin bond layer and good thermal contact. Epoxy will cure in 12 h at 25 °C or in 2 h at 66 °C.

**NOTE:** If using an electrically conductive adhesive or solder, it is important that the excess does not 'creep-up' the edges of the sensor or come in contact with the sensor leads. There is a thin braze joint around the sides of the SD package that is electrically connected to the sensing element. Contact to the sides with any electrically conductive material will cause a short.

3. Follow manufacturers' instructions for adhesive curing schedule. Never heat the sensor above 200 °C.

## LEAD ATTACHMENT

1. Although the Cernox® sensor in a SD package is a two-lead device, measurements should be preferably made using a four-wire configuration to avoid uncertainties associated with lead resistance.  
*Two-lead measurement scheme* — the leads used to measure the voltage are also the current carrying leads. The resultant voltage measured at the instrument is the sum of the temperature sensor voltage and the voltage drop across the two leads. (See Figure 1A.)  
*Four-lead measurement scheme* — the current is confined to one pair of current leads with the sensor voltage measured across the voltage leads. (See Figure 1B.)
2. Lead polarity: None.



**Figure 1. Two-Lead Versus Four-Lead Measurements**

- Strip the insulation from the connecting wires by delicately scraping with a razor blade, fine sand paper, or steel wool. Phosphor-bronze or manganin wire, in sizes 32 or 36 AWG, is commonly used as the connecting lead wire. These wires have low thermal conductivity and high resistivity which help minimize the heat flow through the leads. Typical wire insulation is Polyvinyl Formal (Formvar™) or Polyimide (ML). Formvar™ insulation has better mechanical properties such as abrasion resistance and flexibility. Polyimide insulation has better resistance to chemical solvents and burnout.
- Prepare the connecting wire ends with a RMA (rosin mildly active) soldering flux, tin them with a minimal amount of 60% Sn/40% Pb solder. Use a low wattage soldering iron which will not exceed 200 °C.
- Clean off residual flux with rosin residue remover. The sensor leads can be prepared in an identical manner.
- Join one sensor lead with two of the connector wires. Apply the soldering iron to the connector wire above the joint area until the solders melt, then remove the iron. Repeat for the other set of connector wires and the other sensor lead. Heat sinking the SD package with a flat jaw alligator clip is good practice to eliminate heat build up at the sensor element.
- Avoid putting stress on the device leads and leave enough slack to allow for the thermal contractions that occur during cooling which could fracture a solder joint or lead. Some epoxies and shrink tubing can put enough stress on lead wires to break them.

#### HEAT SINKING/THERMAL ANCHORING

- Since the area being measured is read through the base of the sensor, heat flow through the connecting leads can create an offset between the sensor chip and the true sample temperature. Thermal anchoring of the connecting wires is necessary to assure that the sensor and the leads are at the same temperature as the sample.
- Connecting wires should be thermally anchored at several temperatures between room temperature and cryogenic temperatures to guarantee that heat is not being conducted through the leads to the sensing element. Two different size copper bobbins are available from Lake Shore for heat sinking connecting leads: P/N 9007-900 and 9007-901.
- If connecting wires have a thin insulation such as Formvar™ or Polyimide, a simple thermal anchor can be made by winding the wires around a copper post, bobbin, or other thermal mass. A minimum of five wraps around the thermal mass should provide sufficient thermal anchoring. However, if space permits, additional wraps are recommended for good measure. To maintain good electrical isolation over many thermal cycles, it is good practice to first varnish a single layer of cigarette paper to the anchored area then wrap the wire around the paper and bond in place with a thin layer of IMI 7031 varnish. Formvar wiring insulation has a tendency to craze with the application of IMI varnish. If used, the wires cannot be disturbed until the varnish is fully cured and all solvents have evaporated (typically 12 to 24 h).
- A final thermal anchor at the sample itself is good practice to ensure thermal equilibrium between the sample and temperature sensor.

**CRYOGENIC ACCESSORIES** – Recommended for proper installation and use of CX-10XX-SD package:

**Stycast® Epoxy 2850FT** (P/N 9003-020, 9003-021): Permanent attachment, excellent low temperature properties, poor electrical conductor, low cure shrinkage.

**Apiezon® N Grease** (P/N 9004-020): Low viscosity, easy to use, solidifies at cryogenic temperatures, excellent lubricant.

**IMI 7031 Varnish** (P/N 9009-002): Nonpermanent attachment, excellent thermal conductor, easy to apply and remove.

**Indium Solder** (P/N 9007-002-05): 99.99% pure, excellent electroplating material, foil form.

**Phosphor-Bronze Wire** (P/N 9001-00X): Available in single, duo, and quad strands, no magnetic attraction, low thermal conduction.

**Manganin Wire** (P/N 9001-00X): Low thermal conductivity, high resistivity, no magnetic attraction.

**Heat Sink Bobbin** (P/N 9007-900 Large, 9007-901 Small): Gold-plated oxygen-free high-conductivity (OFHC) copper bobbins.

**Instruments:** Lake Shore sells a complete line of instrumentation used with the Cernox® sensors, such as: current sources, cryopump monitors, temperature controllers, monitors and thermometers, temperature scanners and transmitters.

For complete product description and detailed specifications on the above accessories and instruments, consult the Lake Shore Temperature Measurement and Control Catalog, call (614) 891-2243, e-mail [sales@lakeshore.com](mailto:sales@lakeshore.com), or visit our website [www.lakeshore.com](http://www.lakeshore.com).